





Clesta LED Dental Light

A blend of LED capabilities and ease of use delivers a bright and comfortable work environment.

Easy on the Eyes

A stepless light adjustment provides operators with the ability to fine-tune illumination to fit individual operatory lighting conditions. Adjustment on light intensity creates a comfortable work setting by minimizing eye fatigue.

Positioning Flexibility

A three-axis light head rotation provides operators with the flexibility to direct Illumination precisely where it is needed within the patient's oral cavity.



Color Rendering

The 5 LED bulb array provides cool white light that is excellent for tooth shade matching and differentiation of soft tissues.

Asepsis Benefit

Touchless sensor facilitates hands-free activation for minimizing the risk of cross-contamination of the light head surfaces.





Direct Illumination

Direct LED illumination maximizes efficiency and does not require a cooling fan.



Long LED life and low power consumption adds to operatory efficiency.

Useful Through Obstruction

Shadowless illumination technology allows for uninterrupted lighting to the patient oral cavity, even with a 50% obstruction between the light head and the patient.







Mounting Options

Configurations include ceiling mount, track mount, post mount, curved post unit mount, swing mount.




Belmont Equipment


AL-320 Unit Mount

AL-301 Post Mount

160°

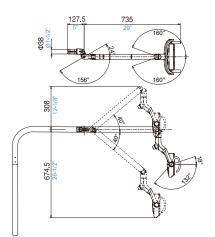

AL-302 Ceiling Mount

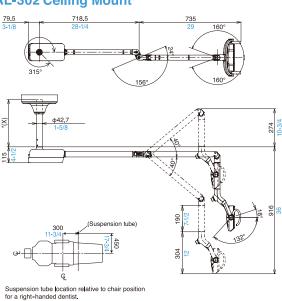


*Select suspension tube length based upon ceiling height

SPECIFICATIONS

Light Source:	Five LEDs
Focal Distance:	25-5/8″ (650 mm)
Intensity:	288 - 2,601 FTC (Stepless)
Composite Safe Mode:	288 FTC
Color Temperature:	5,000° K at any intensity
Color Rendering Index:	85
Light Pattern:	3-3/8″ x 6-1/8″ (85 x 155 mm)
Operation:	Touchless sensor
Intensity Adjustment:	Stepless rotating knob
Head Positioning:	Rotation of head at three axis

Belmont Equipment 101 Belmont Drive Somerset, NJ 08873 (800) 223-1192 www.belmontequip.com Takara Company Canada Ltd. 2076 South Sheridan Way Mississauga, ON L5J 2M4, Canada (800) 268-5351 dental.takarabelmont.ca



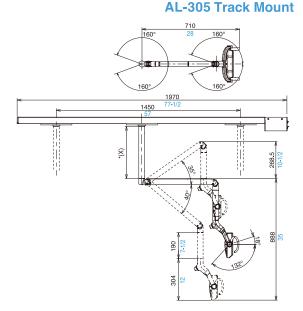

79.5

320.5

718.5

156

Unit : mm