

Forest Dental 7020 PRO Series

Affordable Excellence + Designer Friendly[®]




INTRODUCING THE 7020 PRO SERIES.


**NO MATTER HOW YOU SPIN IT, IT'S ALL
ABOUT REAR-MOUNT SIDE DELIVERY.**

“The ergonomics of the rear-mount side-delivery system are especially important because improper sitting and reaching is what shortens our careers. With this system, you’re sitting in the proper position at 11 o’clock, bent from the hip, and you don’t have to take your eyes away from the patient or reach. Everything is where it needs to be at all times.”

- Brad Boeke, DDS

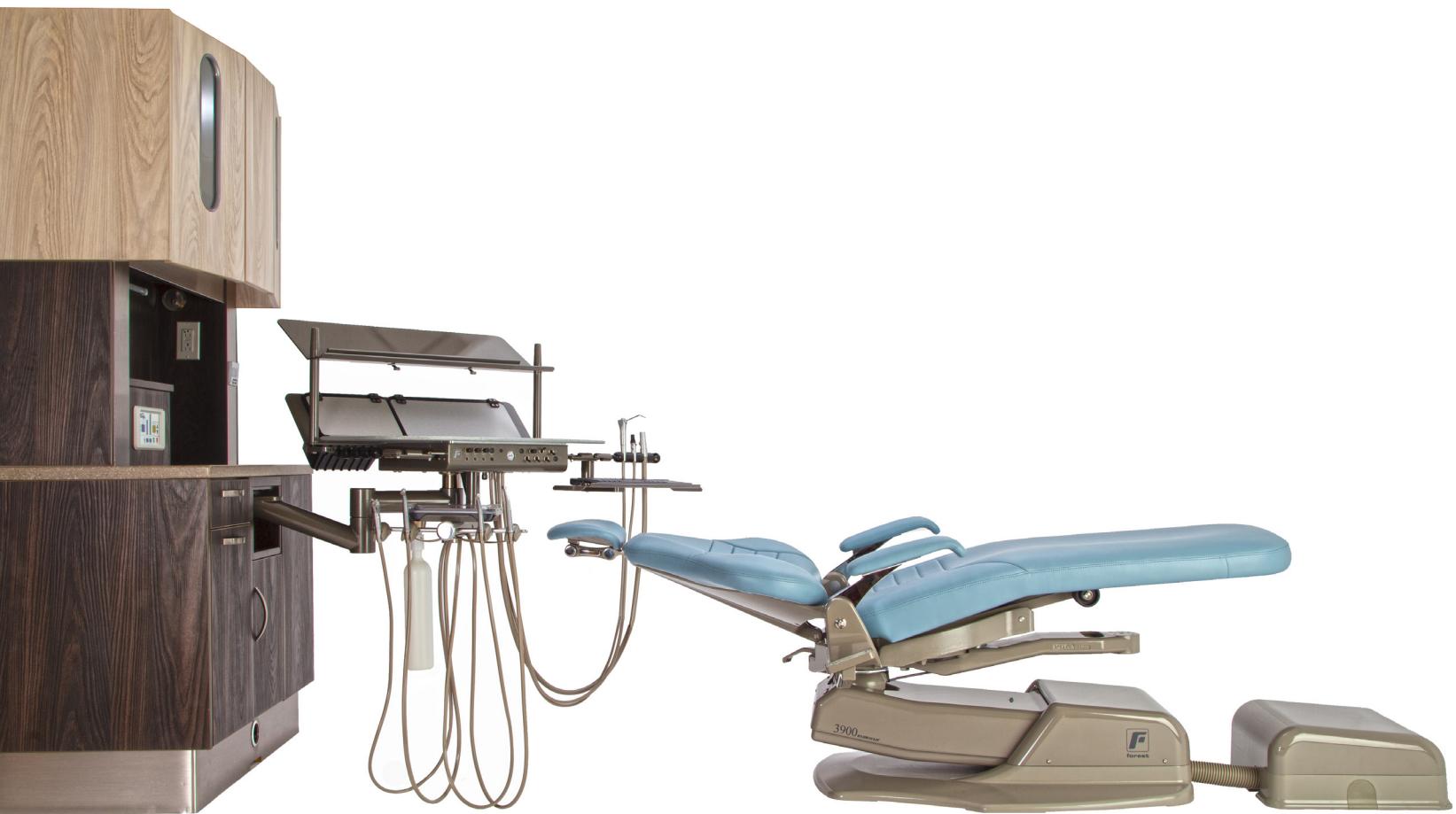



XX

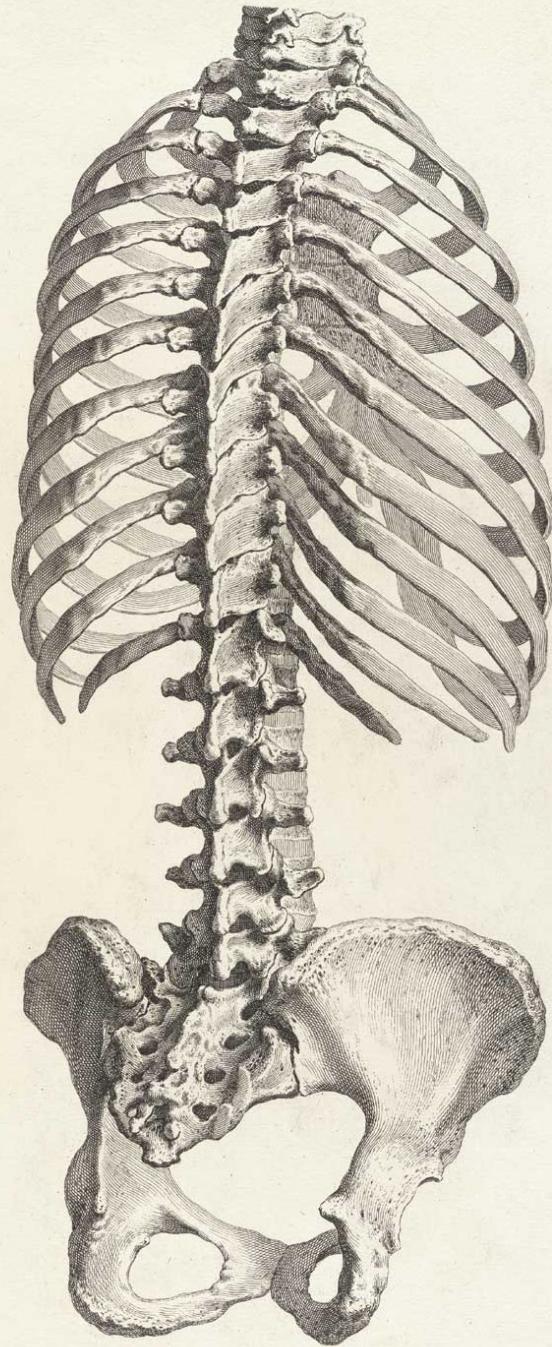

CONSIDER THE SPINE

**DENTISTRY IS DESIGN CENTRIC.
THE 7020 PRO SERIES EMBRACES
ERGONOMICS AND BEAUTY...**

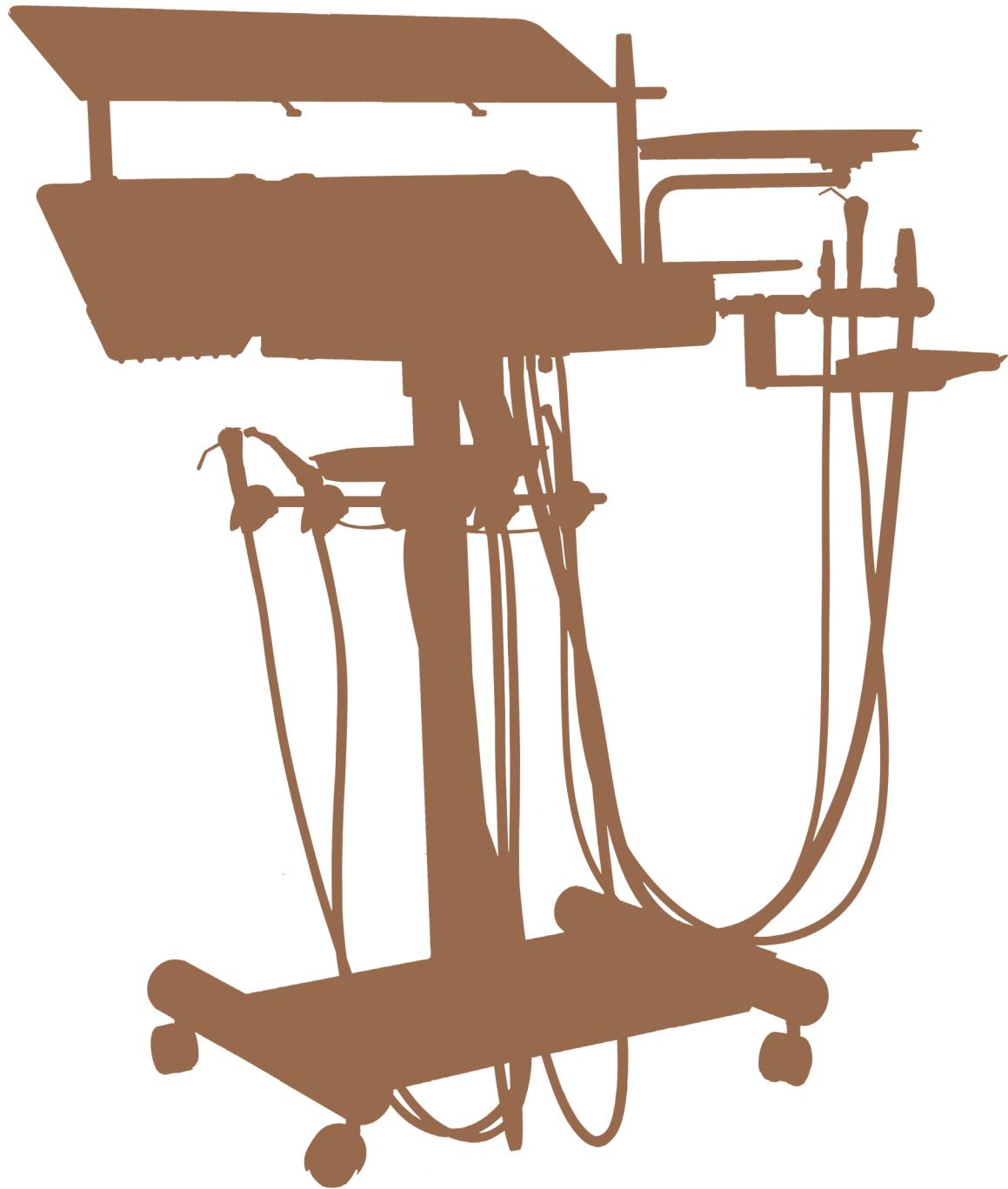

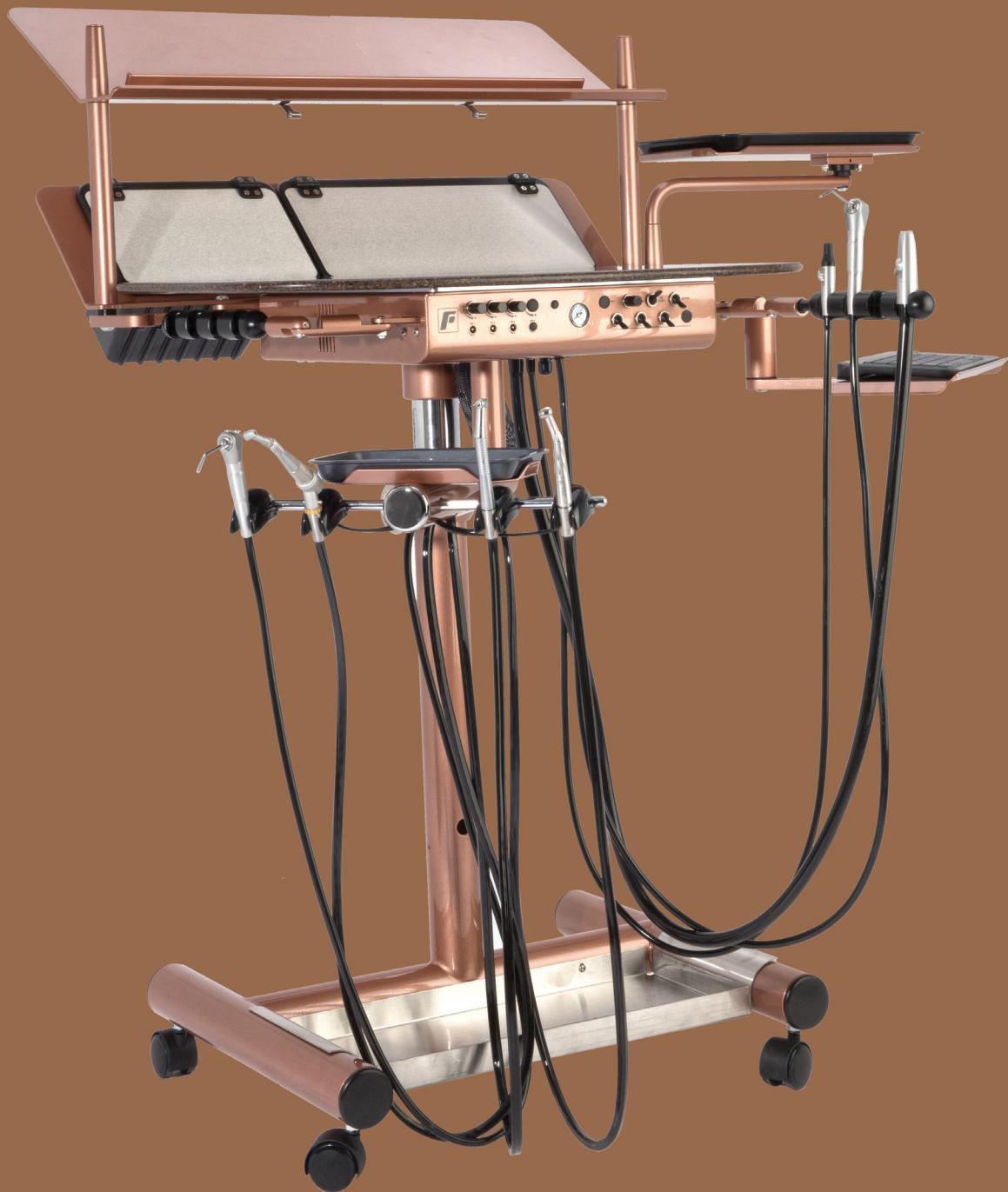



WHY FOREST?

We are very proud to introduce you to our 7020 PRO series side-delivery systems. Compare them with any delivery system—at any price—from any North American manufacturer. You'll easily discover why the 7020 PRO series offers you the highest quality product available today.

The Proof

Comfort, efficiency and durability can be seen—and felt—in every detail done just right. An investment in Forest means:

- More metal, less plastic for a long, stain-free life
- Aluminum chassis and all-metal parts under the hood with life-time warranty
- Four mounting options to best fit your treatment room layout

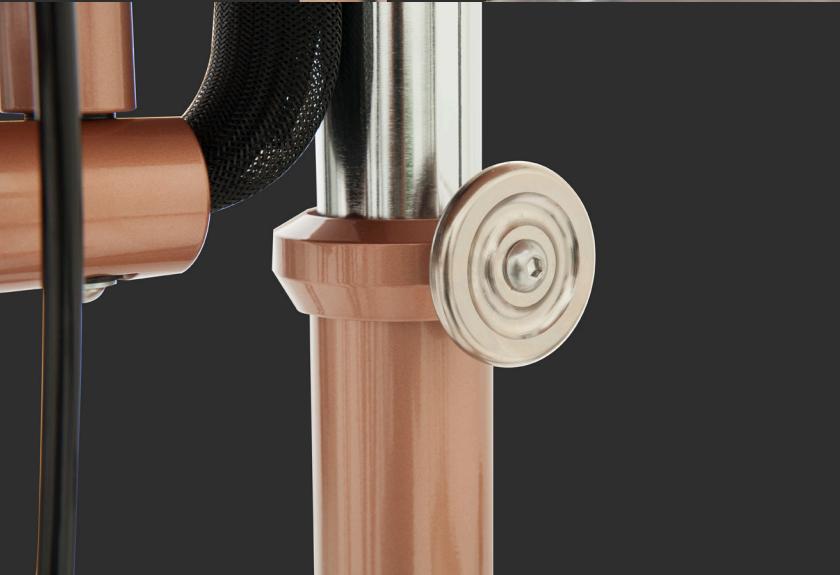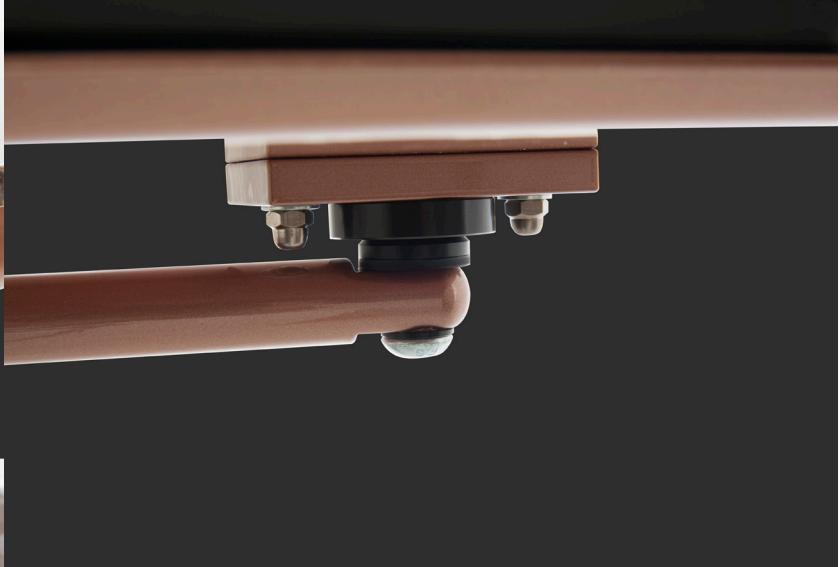

**EVERY DETAIL
DONE RIGHT.**



W E K E E P I T T I G H T

7020 PRO SERIES TOP DOWN FEATURES

Many little things done well. It's all in the details more metal and less plastic for a long beautiful life.

1. Upper Shelf

- Procedure Shelf in 3/16" rust-free aluminum 10 x 32"
- Close to work surface for less reach
- Aluminum radius-ed 3/16" rail support tubs
- Secondary work surface at bottom of shelf

2. Shelf Tower

- 1" round aluminum stock
- Optional tray holders for additional work surface
- Optional monitor mount

3. Bin Rack & Bins

- 3/16" aluminum plate bin rack
- ABS molded plastic bins
- Customizable dividers and shelves
- Brush-anodized aluminum bin lids with quick-release stainless steel spring latch

4. Work Surface

- Choice of Corian® solid surfaces
- Large work surface, finished size 480sq." (16.5" X 30")
- Large aluminum reinforcing plate
- Balanced with 60lb air-spring
- Rounded corners for easy pass-by

5. Chassis

- Life-time warranty
- Aluminum construction
- Large capacity for ancillary options (full 492 sq inches of usable space)
- Powder coat paint or electroless nickel plating on all surfaces
- 4 TIG-welded gusseted corners

6. Key-board tray

- Movable left to right
- Stores under control head
- Structural billet aluminum arm
- Aluminum construction

7. Control

- Life-time warranty on all metal parts
- No water in master control block for long trouble-free life
- Standard hp lock-outs accommodating 4-handed dentistry (allows burr removal with another hp in use)
- Very fine water coolant adjustment
- Tool kit for functionality adjustments
- Room for ancillaries including electric hp, Piezo & magneto strictive scalar, hp light source

8. Burr-Buddy

- Extra work surface for mirror, explorer

9. Chair Control Touchpad (Optional)

- Mounted in middle of hp holder bar



DRESSED UP


DRESSED DOWN



MAKE IT SIMPLE




MOUNT IT

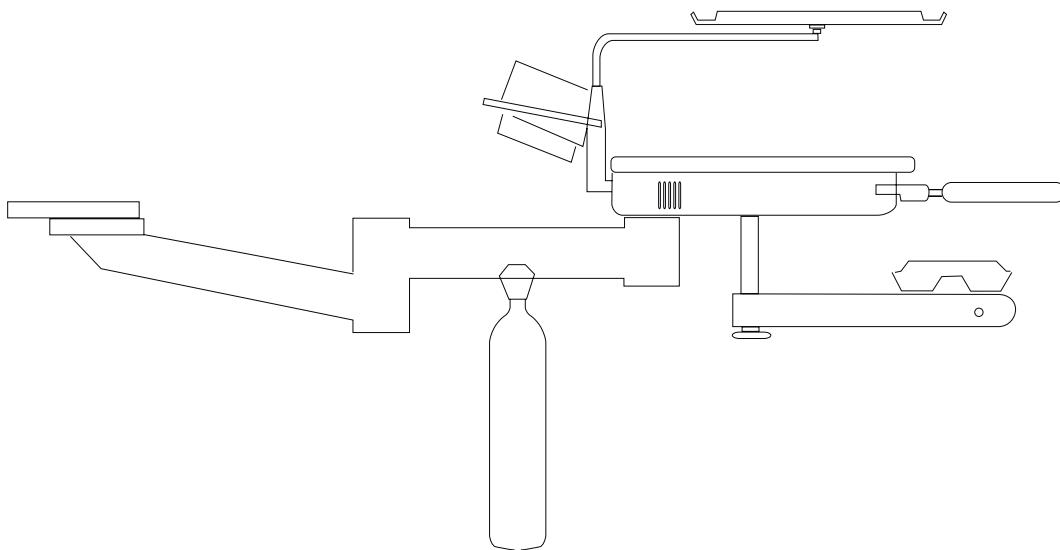


Beneath Counter Mount

The very popular beneath counter mount is anchored by a half-inch steel plate that can be integrated into your custom-made rear cabinet, or as part of a Modular & Custom Cabinetry rear cabinet (PFR8144), which includes the steel plate pre-installed.

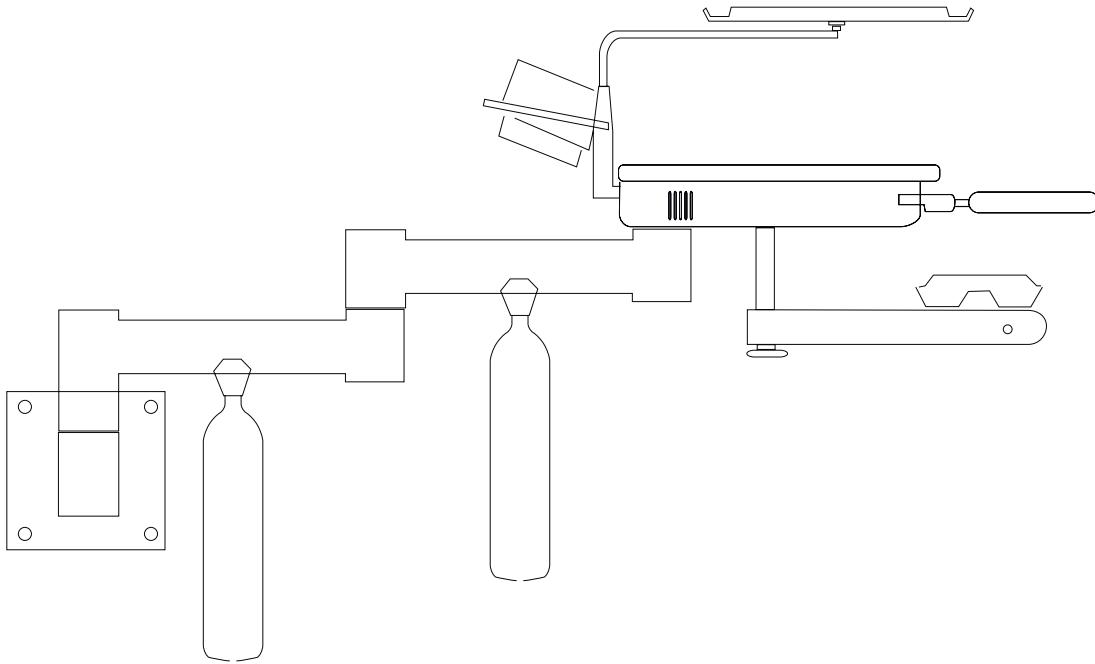

Wall Mount

The wall mount is often used with an upper cabinet module for the storage of tubs and cassettes. The 7020 can replace the lower cabinet because it contains the industry's largest storage space "under the hood" for ancillaries that would otherwise be installed in cabinetry. The money saved by eliminating the lower cabinet goes a long way toward purchasing other practice-enhancing items and is an excellent solution for shorter treatment rooms.

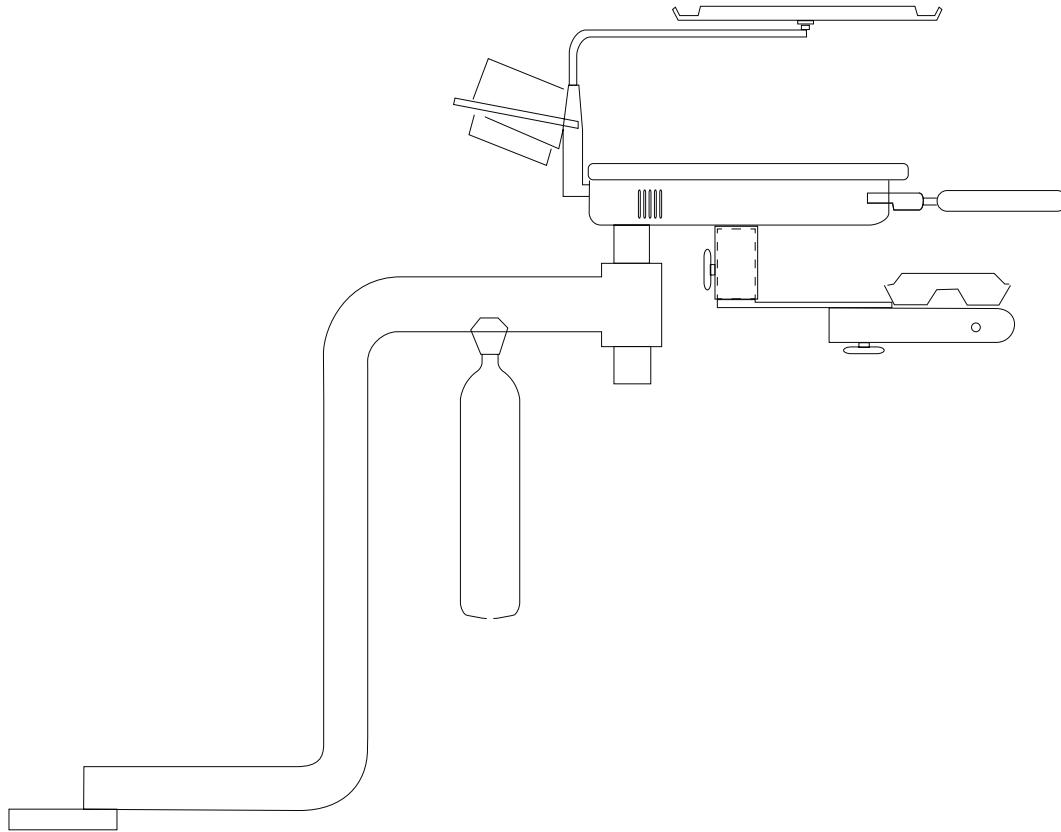


Under Cabinet Mount

The under cabinet mount will free up space for a drawer or nitrous mount by moving the mount underneath the cabinet. Designed to be compatible with all big-brand rear cabinets, the 7020 under cabinet mount is an excellent alternative to their delivery systems.

OR PUT IT ON A CART











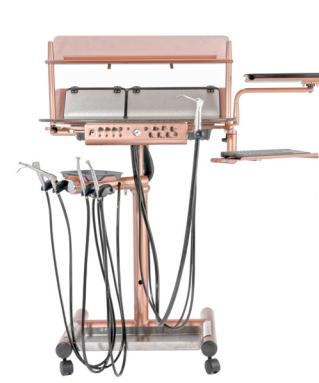


Cart Mount

The cart mount can be ordered with your desired length of umbilical, allowing for a greater range of movement left to right.



TIME TO PLAY

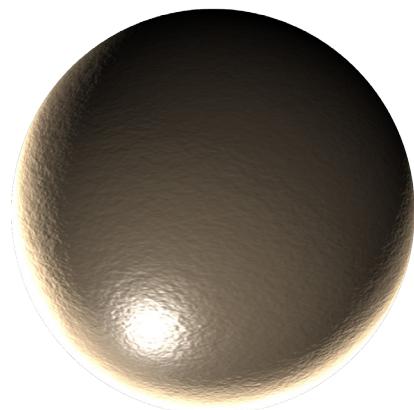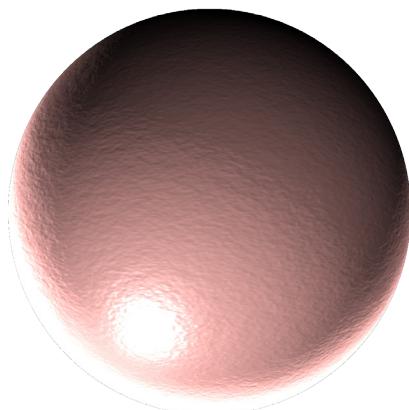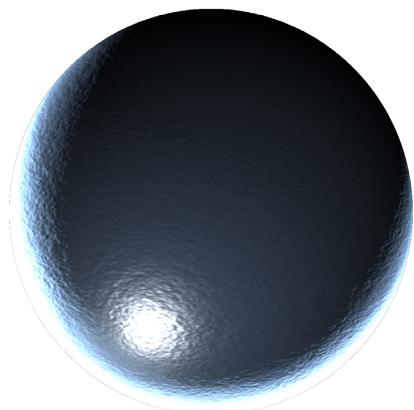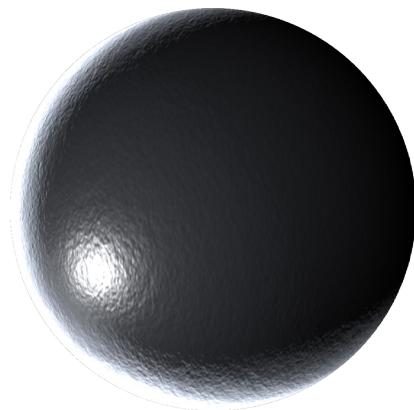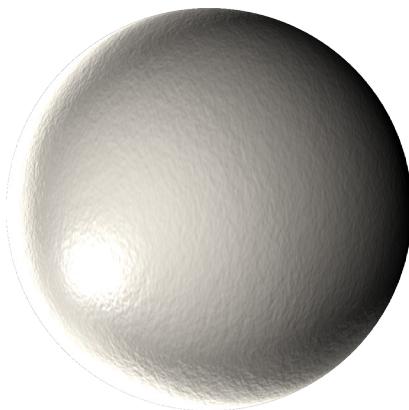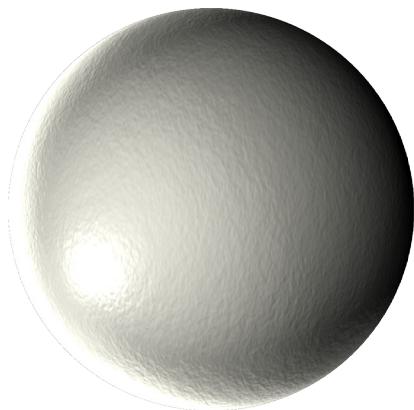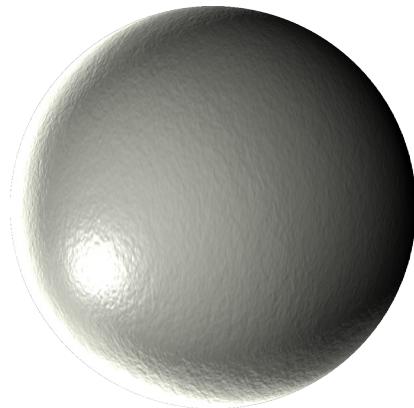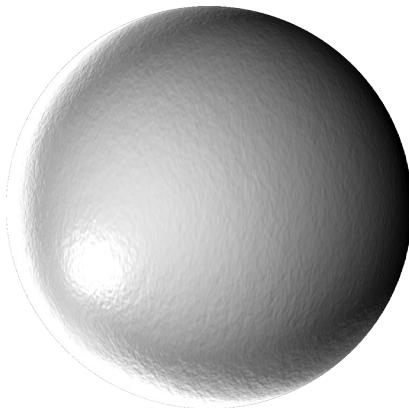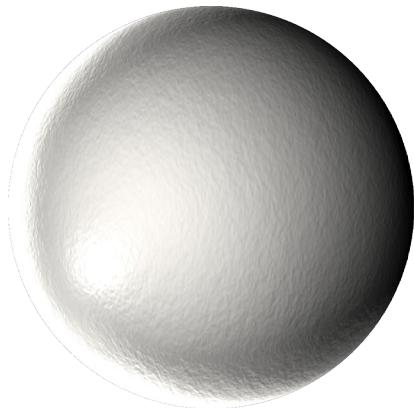



Forest Dental

Phone

- 800 . 423 . 3555
- 503 . 640 . 3012

Fax

- 503 . 693 . 9715

Digital

- sales@forestdental.com
- forestdental.com

Address

- 6200 ne cherry drive hillsboro oregon 97124