



POWERVAC® G DRY VACUUM

The Midmark PowerVac® G provides on-demand energy management that creates a new standard for dental vacuums. Other vacuums operate at full force regardless of demand. The PowerVac G automatically increases the pump speed to match demand, giving you exactly the vacuum you need, no more and no less which reduces energy consumption by as much as 83%.

FEATURES

- Same features as the PowerVac® with additional benefits
- Oil-less and water-less design
- Virtually maintenance free
- Quiet operation
- Compact design uses similar floor space as conventional wet-ring vacuums
- Consistent high performance – up to 18" Hg
- Can replace many types of existing vacuum technologies regardless of plumbing configuration
- Upgradable flexibility – can accommodate 1–20 users
- Unmatched product warranty – 5-year/10,000-hour full-product; 10-year/20,000-hour pump warranty

USES UP TO
83%
LESS ELECTRICITY
THAN COMPETITIVE
SYSTEMS

REDUCES UP TO
30%
IN HEAT EMISSIONS
OVERALL

33%
QUIETER
THAN SIMILIARLY
SIZED DRY VACUUMS

5-YEAR
PRODUCT AND
10-YEAR PUMP
WARRANTY

THE NEXT GENERATION IN DENTAL VACUUM TECHNOLOGY



The powerful 3-HP motor can meet the vacuum demands of up to 10 users simultaneously.





The PowerVac G control panel allows you to easily select the vacuum level that meets your needs.



The PowerVac G has a large-volume separator to deliver the performance demanded by most practices.

POWERVAC® G COMBINES INDUSTRY-LEADING PERFORMANCE WITH AN ENVIRONMENTALLY FRIENDLY, ENERGY-EFFICIENT AND WATER-FREE SYSTEM

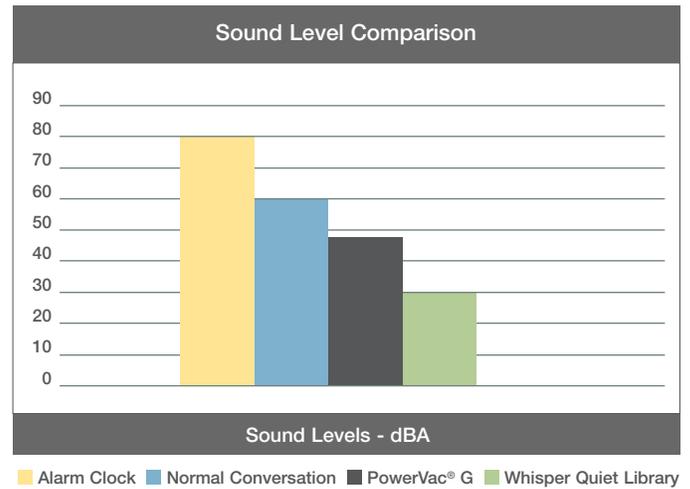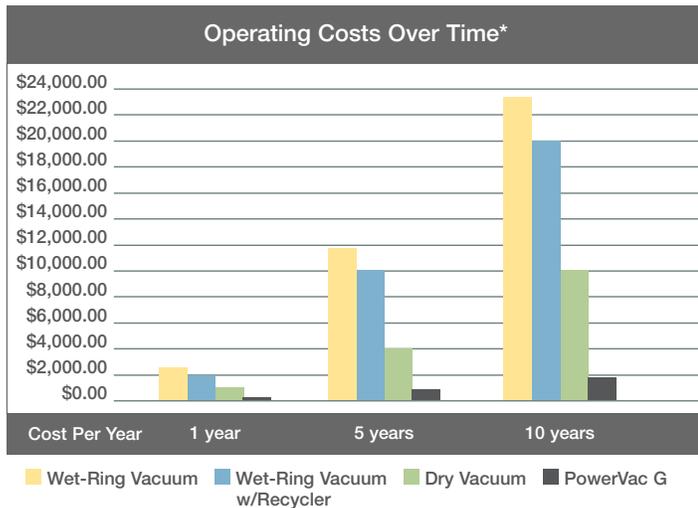

MIDMARK POWERVAC G SPECIFICATIONS

Model/Order No.	G3 (SINGLE)	G5 (SINGLE)	G6 (TWIN)	G7 (SINGLE)	G10 (TWIN)	G14 (TWIN)
No. Users	3-5	5-7	6-10	7-10	10-14	14-20
No. Vacuums	1	1	2	1	2	2
Total Horsepower	3	3	6	3	6	6
Vacuum (" Hg)	up to 18	up to 18	up to 18	up to 18	up to 18	up to 18
Voltage	208-230 V 60 Hz	208-230 V 60 Hz	208-230 V 60 Hz	208-230 V 60 Hz	208-230 V 60 Hz	208-230 V 60 Hz
Total Amps	20	20	20 x 2	20	20 x 2	20 x 2
Circuit Breaker Amps	30	30	30 x 2	30	30 x 2	30 x 2
Dimensions Overall: H x W x D	59" x 22" x 22"	59" x 22" x 22"	59" x 56" x 22"	59" x 22" x 22"	59" x 56" x 22"	59" x 56" x 22"
Base: H x W x D	25" x 22" x 20"	25" x 22" x 20"	25" x 22" x 20"	25" x 22" x 20"	25" x 22" x 20"	25" x 22" x 20"
Tank: H x W x D	34" x 22" x 22"	34" x 22" x 22"	34" x 22" x 22"	34" x 22" x 22"	34" x 22" x 22"	34" x 22" x 22"
Shipping Weight	Base - 431 lb Tank - 90 lb	Base - 431 lb Tank - 90 lb	Base - 431 lb x 2 Tank - 90 lb x 2	Base - 431 lb Tank - 90 lb	Base - 431 lb x 2 Tank - 90 lb x 2	Base - 431 lb x 2 Tank - 90 lb x 2



ACCESSORIES

- 01| Precision Clense Plus Evacuation System Cleaner** keeps your vacuum thoroughly clean, naturally disinfected and deodorized.
- 02| The Solmetex NXT Hg5® amalgam separator** is ISO 11143:2008 certified at greater than 95% separation and meets all "EPA Dental Rule" requirements.

Midmark is an ISO 13485 and ISO 9001 Certified Company. For more information or a demonstration, contact your Midmark dealer or call: 1-800-MIDMARK Fax: 1-877-725-6495 Outside the U.S.A. call: 1-937-526-3662 Fax: 1-937-526-8214 or visit our website at midmark.com

© 2018 Midmark Corporation, Dayton, OH. Products subject to improvement changes without notice. Litho in U.S.A. 007-10037-00 Rev. A1 (2/18)



Complies with NFPA 99C level 3



Because we care.