



INTRAORAL X-RAY

PREVA INTRAORAL DC X-RAY SYSTEM adds superior technology to the most common X-ray unit used by dental practitioners, to:

- Identify cavities and visualize tooth roots
- Examine the status of developing teeth
- Monitor the overall health of the bony area surrounding the tooth



PREVA – FLAWLESS IMAGING AND ENHANCED WORKFLOW

Preva packs many powerful features into one attractive, reliable system.

The Preva is the imaging workhorse of the dental office, used to detect cavities or other problem areas. With exceptional image quality and proven mechanics, Preva units are a pleasure to operate with effortless movement, drift-free positioning and minimum patient stress by getting the images needed with just one shot. Preva offers a powerful combination of features designed to provide maximal flexibility, razor-sharp images and years of productivity. Midmark Preva intraoral X-ray systems use a high-frequency direct current (DC) method of image capture in a range of solutions for your practice needs.

Enjoy the next level of convenience through Preva Plus integration with our ClearVision® sensor. ClearVision® high resolution digital sensors achieve the ultimate in digital convenience, allowing captured images to be viewed with Progeny® Imaging Software or with one of several popular practice management and imaging programs.

FEATURES

- The Preva system is designed for effortless movement and fast one-handed positioning.
- The proprietary dual-arm braking system helps the head stay exactly where you position it without sagging or drifting.
- Preva provides a 0.4 mm focal spot for sharp, well-defined details that allow easy identification of pathology with both film and digital sensors.
- The easy-to-use LCD control panel display with large icons helps to identify target tooth, receptor type and patient size while also offering pre-programmed technique settings.
- The Preva Plus is configured with all the features of the Preva DC and Preva 2.0, plus integration with our ClearVision digital sensor, and built-in storage to keep sensors from being damaged or misplaced.
- All Preva models are available in multiple mounting configurations to best suit your operatory needs: wall, cabinet or mobile.
- Designed and packaged for fast and easy installation.
- Auto-calibration during servicing ensures consistent image quality.

ADDITIONAL FEATURES



For consistently superior diagnostic images, selectable options allow you to determine the techniques you prefer to control according to clinical or procedural needs.

- kV – 60, 65 and 70
- mA – 4, 5, 6 and 7
- Time – 0.02 seconds to 2 seconds



Grow with the Preva 2.0 that gives you a high-speed USB connector so you can go digital when your practice is ready.

PREVA SPECIFICATIONS

Line voltage	110-230 VAC +/- 10% 50 or 60 Hz
Tube potential	60, 65 or 70 kVp
Tube current	4, 5, 6, or 7 mA
Exposure time	20 ms to 2 seconds
Cone length	20 cm (8")
Focal spot	0.4 mm (IEC 60336)
Duty factor	1:15
Shipping weight (wall mount)	95 lb
Applicable documents Federal regulations:	21 CFR, UL, IEC, CE
Warranty	2-year, limited
Options	Exposure hand switch 12" cone 8" rectangular cone Choice of 56", 66", 76" or 82" (total reach) Two-stud mounting plate Metal stud-mounting plate
Minimum ceiling requirement	86.5" (219.7 cm)
Cabinet pass through (alternate)	66" Arm: 26.75" (W) x 47.5" (D) (68 cm x 121 cm) 76" Arm: 36.75" (W) x 47.5" (D) (93.5 cm x 121 cm)
Mobile	Base: 26.9" (W) x 36" (L) (68.3 cm x 91.4 cm) Arm reach: 56" (142 cm) Height: 76.9" (195.3 cm)

SENSOR SPECIFICATIONS

Film size equivalent	SIZE 1	SIZE 2
Dimensions	24 x 27 x 5.5 mm	32 x 43 x 5.5 mm
Active area	30 x 20 mm	36 x 26 mm
Sensor cable	90 cm, 3 m	
Computer interface	High-speed USB 2.0	
Computer requirements	Please see midmark.com for operating system and hardware requirements	

Specifications are subject to change without notice.


ORDERING GUIDE

CLEARVISION® INTEGRATED DIGITAL SENSOR SYSTEM

- XSDN-01** ClearVision® Digital Sensor System Size 1 sensor
XSDN-02 ClearVision® Digital Sensor System Size 2 sensor

ORDERING GUIDE

	PREVA DC	PREVA 2.0	PREVA PLUS
SINGLE STUD MOUNT			
76" Reach #1 Sensor #2 Sensor	P7017	DCD7-G1N	DCD7-G1N/S1 DCD7-G1N/S2
66" Reach #1 Sensor #2 Sensor	P7016	DCD6-G1N	DCD6-G1N/S1 DCD6-G1N/S2
56" Reach #1 Sensor #2 Sensor	P7015	DCD5-G1N	DCD5-G1N/S1 DCD5-G1N/S2
DOUBLE STUD MOUNT			
82" Reach #1 Sensor #2 Sensor	P7018-P	DCD8-G2N	DCD8-G2N/S1 DCD8-G2N/S2
76" Reach #1 Sensor #2 Sensor	P7017-P	DCD7-G2N	DCD7-G2N/S1 DCD7-G2N/S2
66" Reach #1 Sensor #2 Sensor	P7016-P	DCD6-G2N	DCD6-G2N/S1 DCD6-G2N/S2
56" Reach #1 Sensor #2 Sensor	P7015-P	DCD5-G2N	DCD5-G2N/S1 DCD5-G2N/S2
MOBILE			
Mobile	P7017GM	DCDM-G0N/H	DCDM-G0N/H/S1 DCDM-G0N/H/S2

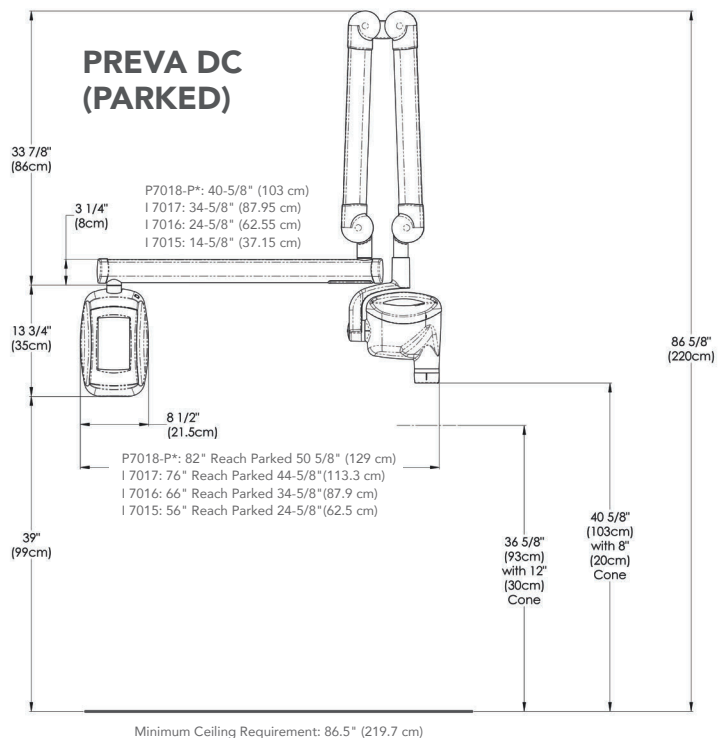

Midmark is an ISO 13485 and ISO 9001 Certified Company.
 For more information or a demonstration, contact your Midmark dealer or call: 1-800-MIDMARK Fax: 1-877-725-6495
 Outside the U.S.A. call: 1-937-526-3662 Fax: 1-937-526-8214 or visit our website at midmark.com

© 2018 Midmark Corporation, Dayton, OH.
 Products subject to improvement changes without notice
 Litho in U.S.A. 007-10044-00 Rev. A1 (2/18)

For additional spec drawings, please review the Preva Dental X-ray System Installation and Service Manual at midmark.com



Because we care.