

# RAMVAC<sup>®</sup> Vacuums





## Bulldog® QT


Rest easy knowing that the Bulldog QT with a 10-year pump limited warranty, will run consistently and efficiently for years to come. Select from two models supporting up to 4 users to meet your specific requirements.



## Bison®

For the power you need in multi-doctor practices, Bison Vacuum Systems are your best choice. Select from four models supporting up to 12 users to meet your specific requirements.

## Rotary Vane Pump Technology

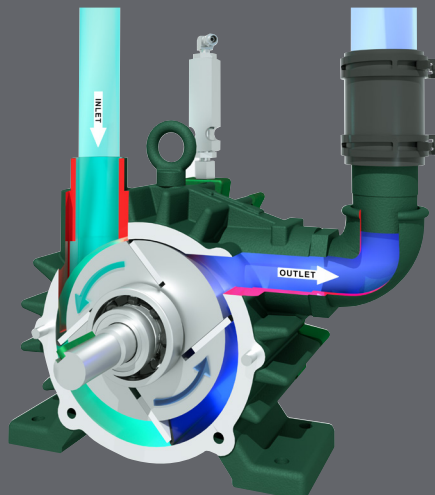

RAMVAC vacuums systems are built with rotary vane pumps. Rotary vane pumps are lubricated and require minimal maintenance but are designed to outlast the service life of non-lubricated systems. The design of a rotary vane pump ensures consistent vacuum performance even as additional HVEs are opened.

## Bulldog QT Combos

Product	Bulldog QT Combo 1	Bulldog QT Combo 2
Part No.	C809001	C809002
Users	3	4
Vacuum Unit	Bulldog QT 1	Bulldog QT 2
Tank Size	15 Gallon (60 L)	15 Gallon (60 L)

## Bison Combos

Product	Bison Combo 3	Bison Combo 5	Bison Combo 7	Bison Combo 9
Part No.	C818003	C818005	C818007	C818009
Users	5	7	10	12
Vacuum Unit	Bison 3	Bison 5	Bison 7	Bison 9
Tank Size	30 Gallon (115 L)	50 Gallon (190 L)	InfiniTank	InfiniTank

## Configurations



Bulldog QT alongside 15 gallon Otter Tank



Bulldog QT stacked on 15 gallon Otter Tank



Bison shown on 30 gallon Otter Tank® using optional utility stand



Bison alongside 50 gallon cylindrical tank

## Vacuum Features

- 100% water free operation
- 10 year pump limited warranty
- Easy and fast annual maintenance\*
- Complete System is UL 2601-1 Listed

\*1 year equals about 2,000 hours of operation. Bulldog QT's S2 controls track cumulative run hours and will display maintenance reminders.



### OWL Touch

S2 controls integrate with the OWL Touch allowing communication between multiple vacuum configurations without costly vacuum equalizers or MPC's.

# Badger™ LF


The Badger LF LubeFree dry vacuum system provides consistent power in the operatory. Worry-free performance and simplicity.



Easily monitor vacuum gauge



Compact design with robust pump

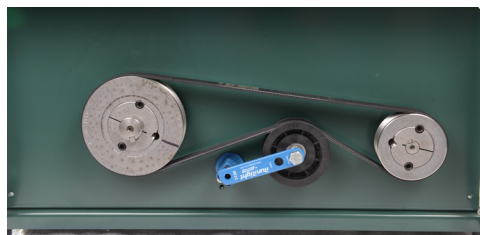

Belt is easily accessible for routine maintenance

## Badger LF Features

- Open frame design means quick maintenance and less down time
- Connects to most pre-existing plumbing, making an easy switch from your old wet-vac system
- Integrated pump, tank and control unit in one compact footprint

## Badger LF Specifications

Product	Badger LF1	Badger LF2
Part No.	BADGLF1	BADGLF2
Users	3	4
Voltage	208-230 VAC, Single Phase	208-230 VAC, Single Phase
Tank Size	5 Gallon	5 Gallon
Full Load	6.7 A	6.7 A
Horsepower	1.5 HP	1.5 HP
Dimension	24.7" W x 18" D x 23.25" H	24.7" W x 18" D x 23.25" H



For more information, please visit [www.DentalEZ.com](http://www.DentalEZ.com)  
Need a RAMVAC part? Visit [www.DentalEZParts.com](http://www.DentalEZParts.com)