


RAMVAC® Utility Room Accessories



Otter Tank®

RAMVAC® tough vacuum tanks, built with flexibility to best fit your needs.

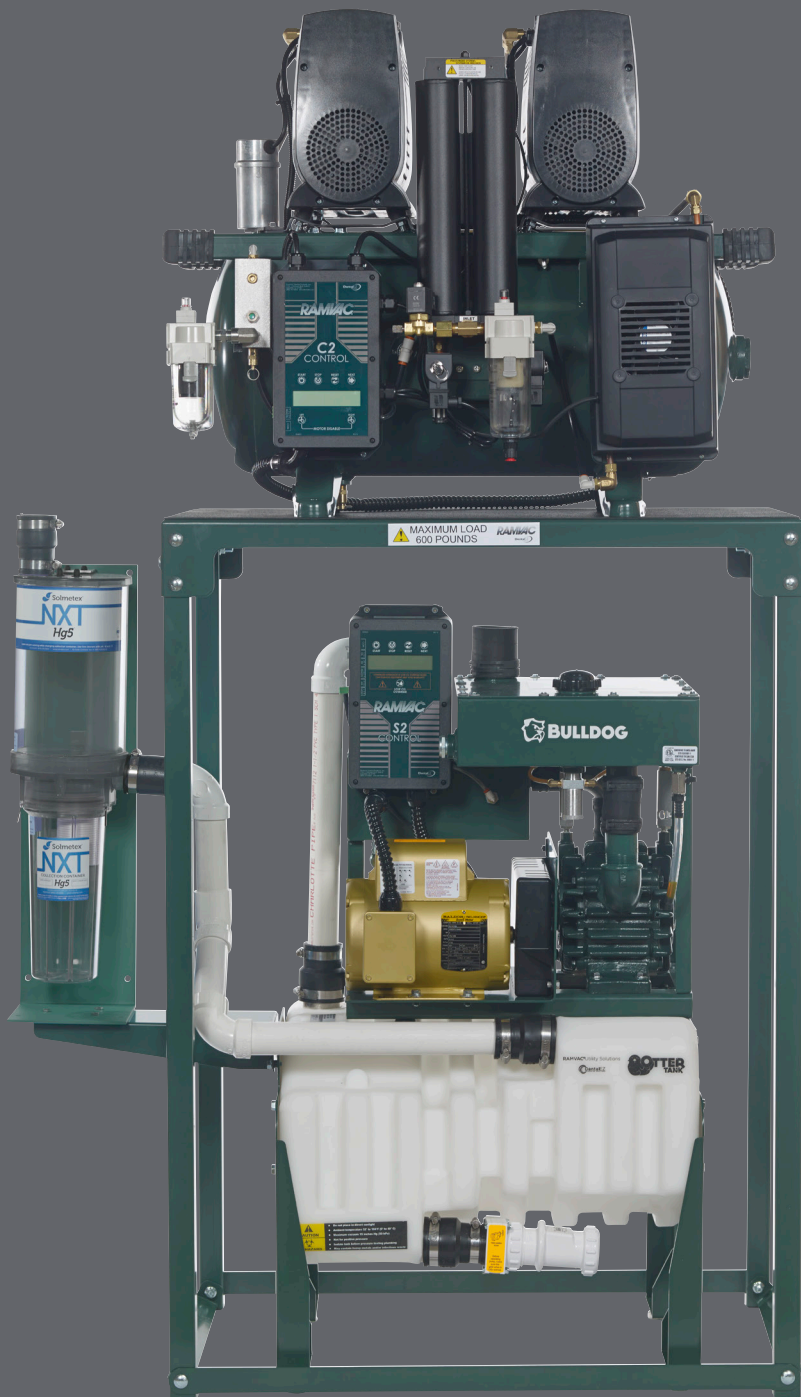

Otter Tank Features

Save Time & Money

- Proven strength and durability
 - Non-corrodible advanced polymer
 - Large, clog resistant drain valve
 - Easy-access top inlet and pump connections
 - Allows multiple “tank+vacuum unit” configurations
 - Easily configured for suspended installation
 - 12” height extender available to raise drain valves
1. Easy access, no glue plumbing connections
 2. Easy access front connection to facility vacuum line
 3. Bulldog mounts directly on top surface. Bison® mounts on optional “Bison Platform”
 4. 12 gauge steel side supports
 5. Cut-outs for routing drain plumbing
 6. Four adjustable non-skid feet
 7. No glue plumbing connections
 8. Proven RAMVAC drain valve
 9. Slide valve
 10. Ribbed construction for added strength
 11. Advanced corrosion-proof lightweight polymer

InfiniTank®

Vacuum separation tanks that automatically drain without the added cost of a pump.



InfiniTank Features

Save Time & Money

- Simple maintenance free design
- Drains without vacuum interruption
- Significantly reduces the size and number of separating tanks
- Completely automated to provide worry free operation
- Intelligent controller records cycle history
- Fits neatly next to walls and in corners

1. Easy access, no glue plumbing connections
2. Low voltage solid-state controller
3. UL60601 listed, medical grade, wall mount transformer, 115V-60Hz/320V-50Hz compatible
4. Ribbed construction and advanced corrosion-proof lightweight polymer for added strength
5. Cut-outs for routing drain plumbing
6. Corrosion proof float switch
7. Proven RAMVAC drain valves
8. Four adjustable feet
9. 12 gauge side supports
10. Reliable, heavy duty, low voltage UL listed solenoid valve for years of trouble free service
11. Bulldog® mounts directly on top surface. Bison® mounts on optional Bison Platform.




OWL Touch

Create a Smart Utility Room

The OWL Touch incorporates a full color touch screen display to provide easy control and monitoring of your dental utility room equipment. Simple, easy to understand screen layouts make navigation a cinch. Automate your utility room operation!

Demo the OWL Touch experience online at: DentalEZ.com/OWL.

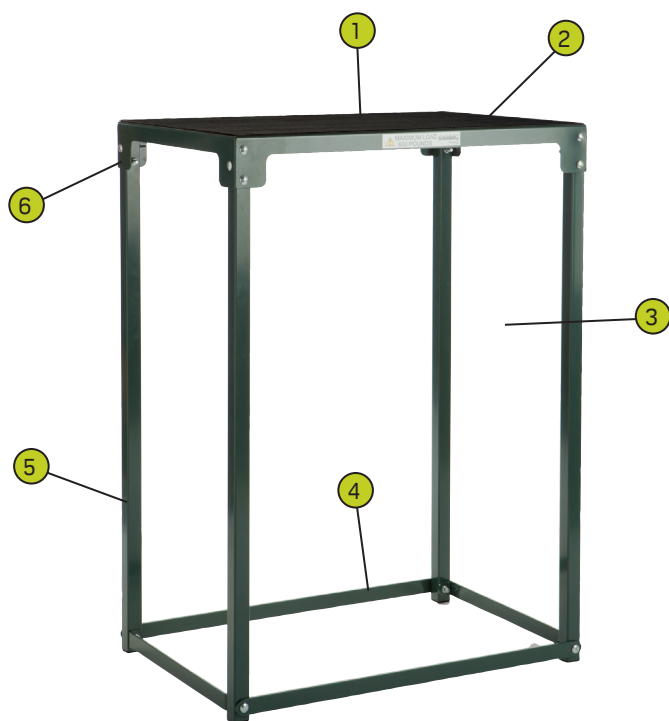

Utility Stands

Maximize Your Space!

1. Load rating is a hefty 600 pounds
2. Rubber mat to reduce noise
3. Wide open front and sides for easy accessibility
4. 12 gauge thick steel straps
5. Sturdy 14 gauge square tubing legs
6. Design provides maximum strength and rigidity



SlugBuster™

Chelators:

- Holds Mercury in suspension for safer discharge to waste
- Softer water for better cleaning

Wetting Agents:

- Specially formulated to allow the cleaner to penetrate faster



For more information, please visit www.DentalEZ.com

Need a RAMVAC part? Visit www.DentalEZParts.com