

BRAVO™

Chamber Autoclaves



THE PERFECT BLEND

of speed and capacity

BRAVO™ is a class B fractionated vacuum sterilizer that incorporates both pre-vac (to remove air) and post-vac (to ensure a dry load every time) cycles, for efficient operation.



Fast

BRAVO™ 17V and 21V are two of the fastest chamber autoclaves on the market, sterilizing and drying instruments in 35-40 minutes.



Gentle

The fresh water feature guarantees fresh water with every cycle. Fresh water reduces the build up of oils and debris.



Versatile

BRAVO™ offers a number of filling & draining options, including a front manual fill and rear automatic fill option that can pull water from a 10L remote bottle placed up to 6 feet away.



Dry

With closed door drying and post vacuum technology, the BRAVO™ is able to dry faster than other chamber autoclaves.



Capacity

BRAVO™ 17V can sterilize 3 full-size cassettes or up to 20 pouches.

BRAVO™ 21V can sterilize 3 oral surgery cassettes or up to 25 pouches.



Cycle Times and Maximum Load Weights

Cycles (134°C @ 4 minutes)	Cycle Times (Including Drying**)		Maximum Load Weight	
	17V	21V	17V	21V
Porous (Handpieces) / Wrapped	35 mins.	40 mins.	*	*
Hollow (Handpieces) / Unwrapped	27 mins.	32 mins.	6.0 kg / 13.6 lbs.	7.5 kg / 16.5 lbs.
Solid / Wrapped	24 mins.	28 mins.	3.0 kg / 6.6 lbs.	4.0 kg / 8.8 lbs.
Solid / Unwrapped	19 mins.	23 mins.	6.0 kg / 13.6 lbs.	7.5 kg / 16.5 lbs.

* 17(V) / 21V: Unpackaged porous materials:1kg /1.25kg, Porous material in single package: 0.75kg / 1.0kg, Solid material/hollow instruments in single pack: 3kg /4kg, Solid material / hollow instruments in double package:1.5kg / 2.0kg.

** Note: Cycle times represent 220V models

Unit Specifications

	BRAVO 17V	BRAVO 21V
Unit size (L x W x H) [†]	22.0 x 19.0 x 16.5" / 560 x 480 x 420mm	25.0 x 19.0 x 16.5" / 635 x 480 x 420mm
Volume	4.5 Gal / 17 L	5.5 Gal / 21 L
Chamber Dimensions (D x L)	10 x 13.5" / 254 x 343mm	10 x 17.5" / 254 x 445mm
Weight	128 lbs. / 55 kg	139 lbs. / 63 kg
Voltage	120V, 20A	220V, 10A

[†] When installed into cabinets, a minimum 4" / 2" (102mm / 51mm) free cool air space is required on the back / sides, and top of the unit.

Recommended Accessories



ENSURE™

Monitoring
for steam sterilization

An easy to use, in-office monitoring
system for steam sterilization.


SCI-BISK-120 - Ensure BI Starter Kit
(25 BI's, Incubator & Log Book)



VistaPure™ AutoFill

Replenish any new or older-model
STATIM and BRAVO autoclave with
distilled-quality water automatically,
enabling truly hands-free convenience.

S3020 - AutoFill for STATIM
S3030 - AutoFill for BRAVO

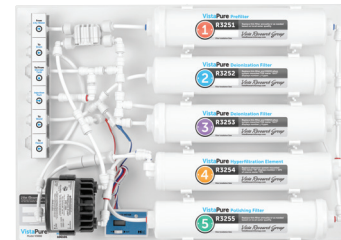

VistaPure™

Water Filtration

Provides two grades of high-quality
water for autoclaves, dental bottles,
ultrasonic cleaners and the rinse
cycle of instrument washers.

V3000

Manufactured by:

SciCan Ltd.

1440 Don Mills Rd., Toronto, ON, Canada, M3B 3P9

Tel: 416-445-1600

Toll free: 1-800-667-7733

customerservice@scican.com

Distributed by:

SciCan Inc.

701 Technology Dr., Canonsburg, PA 15317, USA

Tel: 724-820-1600

Toll Free: 1-800-572-1211

uscustomerservice@scican.com

VistaPure™

are manufactured by:

Vista Research Group

PO Box 321, Ashland, Ohio 44805-0321 USA

Tel: (419) 285-3927

VistaResearchGroup.com